

BAGFILTER PIPET & ROLL

Filter bags with pipet access & built-in closure



The smart
sticker bag!

interscience

Great lab work efficiency!

PATENT PENDING

BagFilter® Pipet & Roll bags, with **BagRack® Slide** storage rack, are an innovative solution for ergonomic work. They are most handy for pathogens research, such as *Listeria* and *Salmonella* analyses, when repeat pipetting from the sample bag is required.

1 Roll

Simply **roll** your bag **3 times** and **close it** with the pre-cut **sticker** for a **waterproof** extra-resistant closure.

2 Slide


Easily **slide your bag** from the side to store your bags during multiple pipettings. Enjoy great **ergonomics** and a **safe, cross-contamination free** work environment.



3 Pipet

Open/close the pipet access easily with the **reclosable sticker**. The filtrate is debris-free thanks to the filter. When multiple pipettings are required for **Pathogens research**, access with pipet tips is made easy by the **lower pipetting height**, without risks of musculoskeletal disorders.

Quick sticker closure

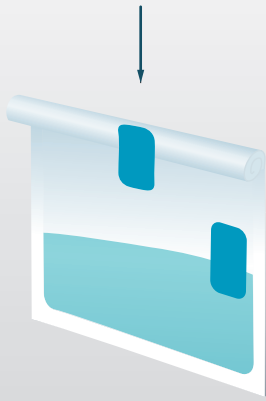

Without any extra tools, simply close with the built-in closure sticker.

BAGFILTER® PIPET & ROLL

Filter bag with built-in closure


With other products

Bag



Closing system

Extra resistant closure

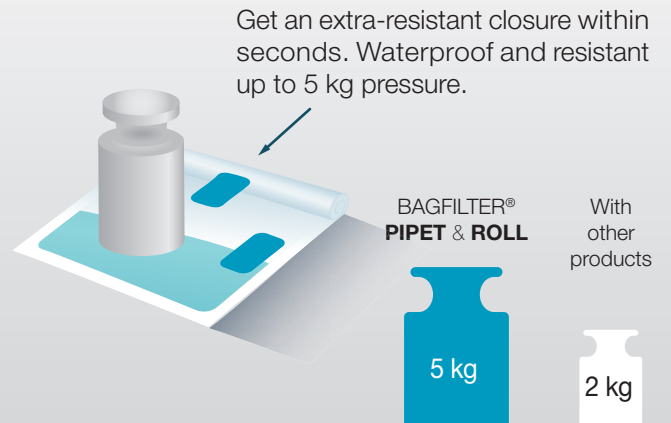

Get an extra-resistant closure within seconds. Waterproof and resistant up to 5 kg pressure.


BAGFILTER® PIPET & ROLL

With other products

5 kg

2 kg

Easy pipetting access

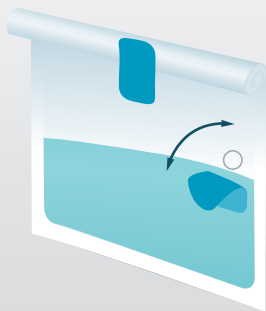

During multiple pipettings for pathogens research, open and close the pipet access easily without opening the bag from the top.

Number of manipulations


BAGFILTER® PIPET & ROLL

1

With other products

1 2 3 4

Slide your bag

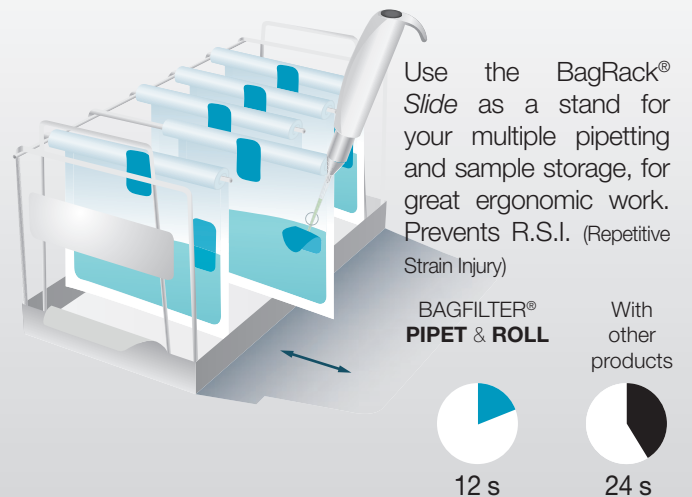

Use the BagRack® Slide as a stand for your multiple pipetting and sample storage, for great ergonomic work. Prevents R.S.I. (Repetitive Strain Injury)


BAGFILTER® PIPET & ROLL

With other products

12 s

24 s

Pathogens research

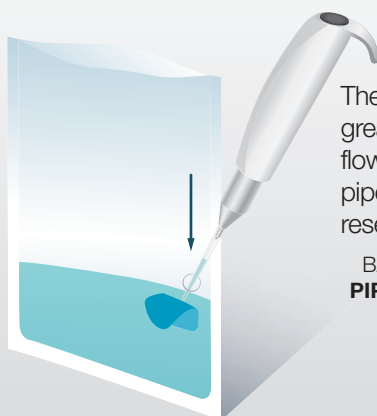

The low pipetting height is great for work under laminar flow cabinets and ideal with pipet tips for Pathogens research.

BAGFILTER® PIPET & ROLL

With other products




10 cm



30 cm

Rapid & hygienic disposal

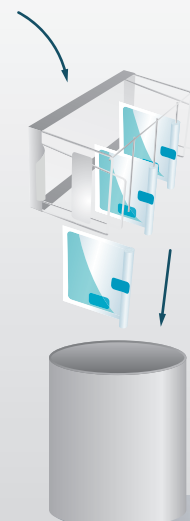

After work, waste disposal is much quicker and without any sample dripping.

Quick

■ ■ ■ ■ BAGFILTER® PIPET & ROLL

■ ■ ■ ■ With other products

Clean

■ ■ ■ ■ BAGFILTER® PIPET & ROLL

■ ■ ■ ■ With other products

3 models for your applications



**BAGFILTER®
PIPET & ROLL**



**BAGFILTER®
PIPET**



**BAGFILTER®
ROLL**

Built-in closure sticker	✓	-	✓
Reclosable easy pipetting sticker	✓	✓	-
Reference	111 720	111 710	111 700
Type of bag	Bag with lateral filter		
Material of the bag	MultiLayer®: reinforced multicoated complex		
Type of filter	Non-woven		
Filter porosity	< 250 microns		
Max volume	400 mL		
Optimal blending volume	50 - 300 mL		
Optimal pipetting volume	250 mL		
Bag dimensions	190 x 300 mm		
Gamma ray treatment	5 to 12 KGy		
In compliance with	ISO 7218 - ISO 6887- FDA-BAM standards		
Pack of	25 bags		
Box of	500 bags		

Ask for free samples

Certified manufacturing

Products made for INTERSCIENCE by Interlab, an ISO 9001 certified company.



Complementary products



DiluFlow®
Gravimetric dilutors
Ref 503 205



BagMixer®
Lab blenders
Ref 025 100



BagRack® Slide **NEW**
Storage rack for bags
Ref 221 440

Your local distributor

interscience

30, chemin du Bois des Arpents - 78860 St Nom - FRANCE
T: +33 (0)1 34 62 62 61 - Email: info@interscience.com

interscience USA & CANADA

32 Cummings Park - Woburn, MA 01801 - USA
P: +1 781 937 0007 - F: +1 781 937 0017 - Email: sales.usa@interscience.com

interscience CHINA

上海市徐汇区肇嘉浜路798号坤阳大厦1903室 - 200030
电话: +86 (0)21-64739390 - +86 189 3097 0733 - 邮址: sales.china@interscience.com

interscience SOUTH-EAST ASIA

The Centropod, 80 Changi Rd - 05-07 Singapore 419715
T: +65 6909 0825 - M: +65 8118 5178 - E-mail: sales.asia@interscience.com

www.interscience.com